
English	INSTRUCTIONS FOR USE: Southern Implants® Drill Extension
Español	INSTRUCCIONES DE USO: Extensión de fresas Southern Implants®
Italiano	ISTRUZIONI PER L'USO: Southern Implants® Prolunga per Punte
Français	MODE D'EMPLOI : Prolongateur de foret Southern Implants®
Deutsch	GEBRAUCHSANWEISUNG: Southern Implants® Bohrer Verlängerung
Português	INSTRUÇÕES DE UTILIZAÇÃO: Southern Implants® Extensão da broca
Ελληνικά	ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ: Southern Implants® Επέκταση τρυπάνου
Svenska	BRUKSANVISNING: Southern Implants® Borrörlängning
Română	INSTRUCȚIUNI DE UTILIZARE: Extensie pentru burghiu Southern Implants®
Türkçe	KULLANMA TALİMATI: Southern Implants® Frez Uzatma

South Africa - Headquarters: 1 Albert Road, Irene, 0062, RSA
T: +27-12-667-1046 | E: info@southernimplants.com

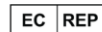Southern Implants Europe AB: Holmgatan 30, S-791 71 Falun, Sweden
T: +46 23 13300 | E: ecrep@southernimplants.com

Subsidiaries

Australia

Southern Implants Australia
T: +61-(0)-8-9466-2627
E: info@southernimplants.com.au

Spain and Portugal

Southern Implants Iberica
T: +34 935 053 507
E: info@southernimplants.es

United Kingdom and Ireland

Southern Implants UK
T: +44-20 8059 4490
E: info@southernimplants.co.uk

USA and Canada

Southern Implants North America Inc.
T: +1-561-472-0990
E: customer care@southernimplants.com

Description

The drill extensions are made of hardened stainless steel. The drill extensions have a latch dimension compatible to ISO 1797. This is in order to connect the drill extension to the handpiece of an implant motor unit. The drill extensions have an internal (where applicable) and external W&H hex connection with dimensions compatible to ISO 1797. This is in order to transmit higher torque values when necessary. The drill extensions are provided non-sterile and are multiple use instruments.

Intended use

The Southern Implants® Drills and Handpiece devices are intended to prepare the implant site/osteotomy and aid in the placement/removal of the endosseous implant and/or prosthetic device(s).

Specifically, the Southern Implants® Drill Extensions are intended to extend the connection of the handpiece device.

Indications for use

The Southern Implants® Drills and Handpiece Devices are indicated for the functional and aesthetic treatment of patients who are eligible for implant-retained prostheses.

Specifically, the Southern Implants® Drill Extensions are indicated in cases where the adjacent structures prevent the drill from reaching the required depth in the bone.

Intended user

Maxillofacial Surgeons, General Dentists, Orthodontists, Periodontists, Prosthodontists, and other appropriately trained and/or experienced medical professionals.

Intended environment

The devices are intended to be used in a clinical environment such as an operating theatre or a dentist consultation room.

Intended patient population

The Southern Implants® Drills and Handpiece devices are intended to be used in patients subject to implant treatment.

Compatibility information

Southern Implants' implants drill extensions are as described in Table A.

Table A –Drill Extensions

Item Code	Material	Description of product
I-DE-G	Stainless Steel (DIN 1.4197)	Drill Extension - Ganter (W&H compatible)
I-DE-GS4	Titanium Grade 5 (ASTM F136), Stainless Steel (17-4PH Condition A)	Drill Extension - Guided Surgery Tool
I-DE-GS5-10	Stainless Steel (DIN 1.4197)	Drill Extension - Guided Surgery Tool, Ø5.2 Sleeve (W&H compatible)
I-DE-K	Stainless Steel (DIN 1.4197)	Drill Extension - Komet
I-DE-MN	Stainless Steel (DIN 1.4305)	Drill Extension - Narrow variation
VXTR	Stainless Steel (DIN 1.4197)	Drill Extension - Versah variation (W&H compatible)

Clinical performance

The Southern Implants® Drill Extensions are designed to assist in the placement or removal of endosseous implants and prosthetic components. The devices are designed to enable secure attachment and controlled torque transfer from the handpiece to the implant or abutment during placement.

Clinical benefits

Southern Implants® Drill Extensions play a vital role in implant-retained prosthetic treatment by enabling secure attachment and controlled torque transfer to the implant or abutment during placement. While their clinical benefits are primarily tied to the overall success of the restorative treatment rather than their direct performance, their role in achieving accurate implant placement supports favourable outcomes. The primary clinical benefits of this treatment approach include the effective restoration of the identified defect and enhanced patient satisfaction. Additionally, successful prosthetic rehabilitation can contribute to improved social interaction, psychosocial well-being, self-esteem, and confidence. Although these benefits are indirect and challenging to quantify, they are closely linked to positive treatment outcomes, such as optimal implant placement, functional restoration, and patient-reported satisfaction with both aesthetics and function.

Storage, cleaning and sterilisation

Limitations on reusable items

A direct value for reusable instruments cannot be given. Frequent processing may have minor effects on the instruments. The product life is normally determined by wear and damage during use, thus instruments if properly cared for and inspected after each use, can be reused many times. Maintain a checklist for these instruments recording the number of uses.

Prior to reprocessing the device, it should be thoroughly inspected and tested to determine its suitability for reuse.

NOTE: during use, handle drills and instruments with sterile tweezers to minimise contamination of the instrument tray and risk of damage to sterile surgical gloves.

Containment

As soon as practically possible, remove all visible residue after use (bone, blood or tissue), by immersing the instrument in cold water (dried soil may be difficult to remove).

Precleaning

Disassemble instruments from hand pieces and all connecting parts from instruments in order to clean soil from obstructed areas. Remove the PEEK bits from placement tools. Rinse with lukewarm water for 3 minutes and remove hardened debris with a soft nylon brush. Avoid mechanical damage during cleaning.

Manual cleaning or automated cleaning

Prepare an ultrasonic bath with suitable detergent (i.e. Steritech instrument cleaner - 5% dilution), sonicate for 20 minutes (alternative methods can be used if proven by the end user). Rinse with purified/sterile water.

NOTE: always follow the instructions for use of the manufacturers for cleaning agents and disinfectants.

Load devices into a thermidisinfector. Run the cleaning and disinfection cycle, followed by the drying cycle.

Drying

Dry the instruments both inside and outside with filtered compressed air or single use, lint free wipes. Pack the instruments as quickly as possible after removal into the storage container. If additional drying is necessary, dry in a clean location.

Inspection

Do a visual inspection of the items to check for damages on instruments.

Packaging

Use the correct packaging material as indicated for steam sterilisation to ensure sterility is maintained. Double packaging is recommended. Where appropriate, the cleaned, disinfected and checked devices can be assembled and placed in instrument trays as applicable. Instrument trays can be double wrapped or placed in sterilisation bags.

Sterilisation

Southern Implants® recommends one of the following procedures to sterilise the device prior to use:

1. prevacuum sterilisation method: steam sterilise the device at 132°C (270°F) at 180 - 220 kPa for 4 minutes. Dry for at least 20 minutes in the chamber. Only an approved wrap or pouch for steam sterilisation must be used.
2. for users in the USA: prevacuum sterilisation method: wrapped, steam sterilise at 135°C (275°F) at 180 - 220 kPa for 3 minutes. Dry for 20 minutes in the chamber. Use a wrap or pouch that is cleared for the indicated steam sterilisation cycle.

NOTE: users in the USA must ensure that the steriliser, wrap or pouch, and all steriliser accessories are cleared by the FDA, for the intended sterilisation cycle.

Storage

Maintain packaging integrity to ensure sterility in storage. Packaging should be completely dry before storage to avoid corrosion and degradation of cutting edges.

Contraindications

The contraindications of all device groups used as part of the specific treatment or procedure apply. Therefore, the contraindications of the systems/medical devices utilized as part of implant surgery/therapy should be noted and the relevant documents consulted.

Contraindications specific to the Southern Implants® Drills and Handpiece Devices:

- Allergies or hypersensitivity to chemical ingredients in the following materials: Titanium, Aluminium, Vanadium, Stainless Steel.

Warnings and precautions

IMPORTANT NOTICE: THESE INSTRUCTIONS ARE NOT INTENDED AS A SUBSTITUTE FOR ADEQUATE TRAINING

- To ensure the safe and effective use of dental implants, new technologies/systems and the Metal Abutment devices, it is strongly recommended that specialised training be undertaken. This training should include hands-on methods to gain competency on the proper technique, the system's biomechanical requirements and radiographic evaluations required for the specific system.
- Improper technique can result in implant failure, damage to nerves/vessels and/or loss of supporting bone.
- Use of the device with incompatible or non-corresponding devices can result in the poor performance or failure of the device.
- When handling devices intraorally, it is imperative that they are adequately secured to prevent aspiration, as aspiration of products may lead to infection or physical injury.
- The use of non-sterile items can lead to secondary infections of the tissue or transfer infectious diseases.
- Failure to adhere to appropriate cleaning, re-sterilization, and storage procedures as outlined in the Instructions for Use can result in device damage, secondary infections, or patient harm.
- Exceeding the number of recommended uses for reusable devices can result in device damage, secondary infection or patient harm.
- The use of blunt drills may cause damage to the bone, potentially compromising osseointegration.

It is crucial to emphasize that training should be undertaken by both new and experienced implant users prior to utilizing a new system or attempting a new treatment method.

Patient Selection and Preoperative Planning

A comprehensive patient selection process and meticulous preoperative planning are essential for successful implant treatment. This process should involve consultation among a multidisciplinary team, including well-trained surgeons, restorative dentists, and laboratory technicians.

Patient screening should include, at a minimum, a thorough medical and dental history, as well as visual and radiological inspections to assess the presence of adequate bone dimensions, the positioning of anatomical landmarks, the presence of unfavourable occlusal conditions, and the periodontal health status of the patient.

For successful implant treatment, it is important to:

1. Minimize trauma to the host tissue, thereby increasing the potential for successful osseointegration.
2. Accurately identify measurements relative to radiographic data, as failure to do so may lead to complications.
3. Be vigilant in avoiding damage to vital anatomical structures, such as nerves, veins, and arteries. Injury to these structures may result in serious complications, including ocular injury, nerve damage, and excessive bleeding.

The responsibility for proper patient selection, adequate training and experience in implant placement, and the provision of appropriate information required for informed consent rests with the practitioner. By combining thorough screening of prospective implant candidates with a practitioner possessing a high level of competence in the use of the system, the potential for complications and severe side effects can be significantly reduced.

High Risk Patients

Special care should be exercised when treating patients with local or systemic risk factors that may adversely affect the healing of bone and soft tissue or otherwise increase the severity of side effects, the risk of complications, and/or the likelihood of implant failure. Such factors include:

- poor oral hygiene
- history of smoking/vaping/tobacco use
- history of periodontal disease
- history of orofacial radiotherapy**
- bruxism and unfavourable jaw relations
- use of chronic medications that may delay healing or increase the risk of complications including, but not limited to, chronic steroid therapy, anti-coagulant therapy, TNF- α blockers, bisphosphonate and cyclosporin

*** The potential for implant failure and other complications increases when implants are placed in irradiated bone, as radiotherapy can lead to progressive fibrosis of blood vessels and soft tissue (i.e., osteoradionecrosis), resulting in diminished healing capacity. Contributing factors to this increased risk include the timing of implant placement in relation to radiation therapy, the proximity of radiation exposure to the implant site, and the radiation dosage at that site.*

Should the device not operate as intended, it must be reported to the manufacturer of the device. The contact information or the manufacturer of this device to report a change in performance is: sicomplaints@southernimplants.com.

Side effects

The clinical outcome of treatment is influenced by multiple factors. The following side effects and residual risks are associated with this device group or the surgical procedures involved in its use. These may require further treatment, revision surgery, or additional consultations with the relevant medical professional.

- Allergic or hypersensitivity reaction(s)
- Anaesthesia, paraesthesia, hyperesthesia, and hypoesthesia (transient or permanent)
- Anatomical landmark damage
- Bruising
- Dental injury during surgery
- Fat emboli formation
- Fistula formation
- Gingival injury
- Infection (acute and/or chronic)
- Localized inflammation
- Loss or damage to adjacent dentition

- Mild thermal injury to adjacent structures (e.g. burn to the lip)
- Nasal bleeding (Zygomatic Implant procedures)
- Nerve damage
- Pain, discomfort or tenderness
- Perforation of the labial or lingual plates
- Peri-implantitis, peri-mucositis or otherwise poor peri-implant soft tissue health
- Sinus Perforation
- Soft tissue irritation
- Surgical side effects such as pain, inflammation, bruising and mild bleeding
- Thermal injury to bone resulting in insufficient levels of osseointegration
- Wound dehiscence or poor healing

Additionally, the normal side effects associated with anaesthesia should also be expected. Note: these side effects and residual risks can vary in both severity and frequency.

Precaution: maintaining sterility protocol

The devices are supplied clean but not sterile in a peel pouch or blister base with peel-back lid. Labelling information is located on the bottom half of the pouch or on the surface of the peel-back lid.

Notice regarding serious incidents

Any serious incident that has occurred in relation with the device must be reported to the manufacturer of the device and the competent authority in the member state in which the user and/or patient is established.

The contact information for the manufacturer of this device to report a serious incident is as follows:

sicomplaints@southernimplants.com.

Materials

Drill Extensions: Stainless Steel according to DIN 1.4197/DIN 1.4305 or 17-4PH Condition A;
 Titanium alloy Ti-6Al-4V (90% titanium, 6% aluminium, 4% vanadium) according to ASTM F136 and ISO 5832-3

Disposal

Disposal of the device and its packaging: follow local regulations and environmental requirements, taking different contamination levels into account. When disposing of spent items, take care of sharp drills and instruments. Sufficient PPE must be used at all times.

Summary of Safety and Clinical Performance (SSCP)

As required by the European Medical Device Regulation (MDR; EU2017/745), a Summary of Safety and Clinical Performance (SSCP) is available for perusal with regard to Southern Implants® product ranges.

The relevant SSCP can be accessed at <https://ec.europa.eu/tools/eudamed>.

NOTE: the above website will be available upon the launch of the European Database on Medical Devices (EUDAMED).

Disclaimer of liability

This product is part of the Southern Implants® product range and should only be used with the associated original products and according to the recommendations as in the individual product catalogues. The user of this product has to study the development of the Southern Implants® product range and take full responsibility for the correct indications and use of this product. Southern Implants® does not assume liability for damage due to incorrect use. Please note that some Southern Implants® products may not be cleared or released for sale in all markets.


Basic UDI

Product	Basic-UDI Number
Basic-UDI for Drills & Handpiece Devices	6009544038759C

Related literature and catalogues

CAT-2004 - Tri-Nex Implants Product Catalogue
 CAT-2005 - IT Implants Product Catalogue
 CAT-2020 - External Hex Implants Product Catalogue
 CAT-2042 - Deep Conical Implants Product Catalogue
 CAT-2043 - Internal Hex Implants Product Catalogue
 CAT-2060 - PROVATA® Implants Product Catalogue
 CAT-2069 - Deep Conical INVERTA® Implants Product Catalogue
 CAT-2093 - Single Platform (SP1) Implants Product Catalogue
 CAT-2095 - External Hex INVERTA® Implants Product Catalogue
 CAT-2089M - SIREAL Guided Surgery System Product Catalogue

Symbols and warnings

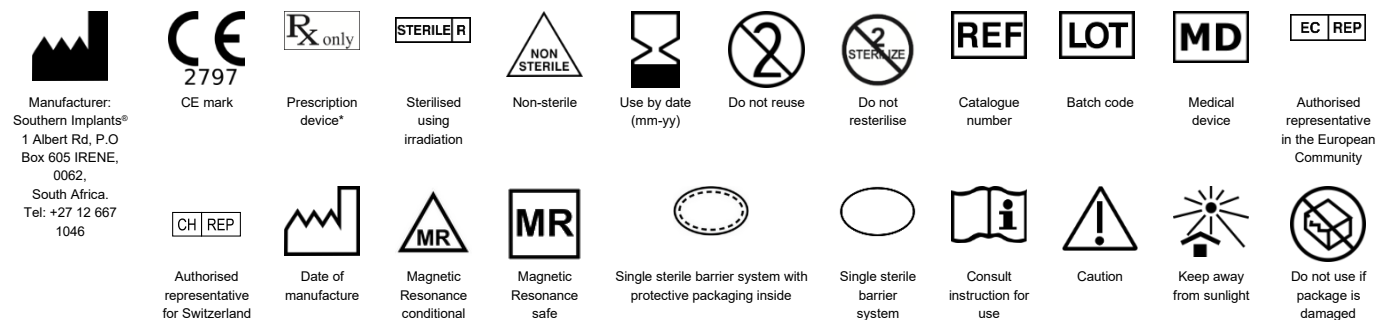* Prescription device: Rx only. Caution: Federal Law restricts this device to sale by or on the order of a licenced physician or dentist.

Canada licence exemption: Please note that not all products may have been licensed in accordance with Canadian law.

All rights reserved. Southern Implants®, the Southern Implants® logotype and all other trademarks used in this document are, if nothing else is stated or is evident from the context in a certain case, trademarks of Southern Implants®. Product images in this document are for illustration purposes only and do not necessarily represent the product accurately to scale. It is the responsibility of the clinician to inspect the symbols that appear on the packaging of the product in use.

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