
English	INSTRUCTIONS FOR USE: Southern Implants® Screwdrivers - Handheld, Wrench insert and Handpiece - Latch
Español	INSTRUCCIONES DE USO: Destornilladores - de mano, inserto para llave y pieza de mano - con pestillo Southern Implants®
Italiano	ISTRUZIONI PER L'USO: Southern Implants® Cacciaviti - Manipolo, inserto per Chiave e Manipolo - Dispositivo di Chiusura
Français	MODE D'EMPLOI : Tournevis Southern Implants® - Portable, insert de clé et pièce à main - Loquet
Deutsch	GEBRAUCHSANWEISUNGEN: Southern Implants® Schraubendreher - Handgerät, Schraubenschlüssel-Einsatz und Handstück - Verriegelung
Português	INSTRUÇÕES DE UTILIZAÇÃO: Southern Implants® Chaves de parafusos - portáteis, de inserção e de peça de mão - engate
Ελληνικά	ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ: Southern Implants® Κατσαβίδια - χειρός, ένθετο κλειδί και χειρολαβή - μάνδαλο
Svenska	ANVISNINGAR FÖR ANVÄNDNING: Southern Implants® Skruvmejslar - handhållen, skiftnyckel inlägg och handstycke - spärr
Română	INSTRUCȚIUNI DE UTILIZARE: Șurubelnițe - portabile, mandrine și partea de mână - dispozitiv de prindere Southern Implants®

South Africa - Headquarters: 1 Albert Road, Irene, 0062, RSA
T: +27-12-667-1046 | E: info@southernimplants.com

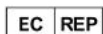Southern Implants Europe AB: Holmgatan 30, S-791 71 Falun, Sweden
T: +46 23 13300 | E: ecrep@southernimplants.com

Subsidiaries

Australia

Southern Implants Australia
T: +61-(0)-8-9466-2627
E: info@southernimplants.com.au

Spain and Portugal

Southern Implants Iberica
T: +34 935 053 507
E: info@southernimplants.es

United Kingdom and Ireland

Southern Implants UK
T: +44-20-8899-6845 / 6 / 7
E: info@southernimplants.co.uk

USA and Canada

Southern Implants North America Inc.
T: +1-561-472-0990
E: customercare@southernimplants.com

Description

The screwdrivers are intended to be used to apply torque to a prosthetic screw or abutment to affix the screw or abutment to the dental implant, or loosen and remove the component from the dental implant or abutment. The screwdrivers have a driving interface on one end that is compatible with the corresponding screw or abutment. This interface allows for pickup of the corresponding component and prevents relative rotation between the screw or abutment and the screwdriver when torque is applied. The screwdrivers are available in various lengths (short 5 mm, medium 14 mm, long 21 mm) to allow for use in both the anterior and posterior where available space may differ.

Intended use

Drills and handpieces are intended to be used to prepare the osteotomy for implant placement.

Indications for use

Dental implant therapy is indicated for patients with one or more missing teeth that wish to have the tooth/teeth restored with a dental implant and prosthetic tooth/teeth.

Intended user

Maxillo-facial Surgeons, General Dentists, Orthodontists, Periodontists, Prosthodontists and other appropriately trained and experienced implant users.

Intended environment

The devices are intended to be used in a clinical environment such as an operating theatre or a dentist consultation room.

Intended patient population

The intended patient population for the device is not dissimilar to that of dental implant therapy. The intended patient population for implant therapy is partially or fully edentulous patients requiring prosthetic dental restoration in the upper or lower jaw. Restorations may comprise single teeth, partial or full bridges and are fixed restorations.

Compatibility information

Southern Implants' implants should be restored with Southern Implants' components. In the Southern Implants' range there are 6 implant connections. Table A indicates the use of Southern Implants screw drivers.

Table A

External Hex						
		Hexed	Unigrip	Quad	Slotted	Torx
	COVER SCREW DRIVERS	HEALING ABUTMENT & PROSTHETIC SCREW DRIVERS	PROSTHETIC SCREW DRIVERS (incl. ASC)			
Handheld	I-CS-HD/L	I-HD-S/M/L I-HD-22U-S/M/L	I-UGI-S/M/L	I-QDI-S/M/L	I-BD-S/M/L	I-HD-ASC-M/L
Handpiece	I-HHD-09	I-HHD-22S/M/L I-HHD-22U-S/M/L	I-HUG-S/M/L	I-HQD-S/M/L	I-HBDS/M/L	I-HASC-M/L
Wrench Insert	I-WI-09	I-WI-22S/M/L I-WI-22U-S/M/L	I-WI-UG-S/M/L	I-WI-QS/M/L	I-WI-BS/M/L	I-WI-ASC-M/L
DC (Deep Conical)						
	Hexed		Torx			
	HEALING ABUTMENT & PROSTHETIC SCREW DRIVERS		ASC PROSTHETIC SCREW DRIVERS			
Handheld	I-HD-S/M/L I-HD-22U-S/M/L		I-HD-ASC-M/L			
Handpiece	I-HHD-22S/M/L I-HHD-22U-S/M/L		I-HASC-M/L			
Wrench Insert	I-WI-22S/M/L I-WI-22U-S/M/L		I-WI-ASC-M/L			
Internal Hex (PROVATA/M-Series)						
	Hexed		Quad		Torx	

	COVER SCREW DRIVERS, HEALING ABUTMENT & PROSTHETIC SCREW DRIVERS	GOLD PROSTHETIC SCREW DRIVERS	ASC PROSTHETIC SCREW DRIVERS
Handheld	I-HD-S/M/L I-HD-22U-S/M/L	I-QDI-S/M/L	I-HD-ASC-M/L
Handpiece	I-HHD-22S/M/L I-HHD-22U-S/M/L	I-HQD-S/M/L	I-HASC-M/L
Wrench Insert	I-WI-22S/M/L I-WI-22U-S/M/L	I-WI-QS/M/L	I-WI-ASC-M/L
TRI-NEX			
	Unigrip	Torx	
	COVER SCREW DRIVERS, HEALING ABUTMENT & PROSTHETIC SCREW DRIVERS	ASC PROSTHETIC SCREW DRIVERS	
Handheld	I-UGI-S/M/L	I-HD-ASC-M/L	
Handpiece	I-HUG-S/M/L	I-HASC-M/L	
Wrench Insert	I-WI-UG-S/M/L	I-WI-ASC-M/L	
IT (Internal Octagon)			
	Torx	Torx	
	COVER SCREW DRIVERS, HEALING ABUTMENT & PROSTHETIC SCREW DRIVERS	ASC PROSTHETIC SCREW DRIVERS	
Handheld	I-SCS-S/M/L	I-HD-ASC-M/L	
Handpiece	I-HSCS-S/M/L	I-HASC-M/L	
Wrench Insert	I-WI-SCS-S/M/L	I-WI-ASC-M/L	
SP1 (Single Platform)			
	Hexed		
	COVER SCREW, HEALING ABUTMENT & PROSTHETIC SCREW DRIVERS		
Handheld	I-HD-S/M/L I-HD-22U-S/M/L		
Handpiece	I-HHD-22S/M/L I-HHD-22U-S/M/L		
Wrench Insert	I-WI-22S/M/L I-WI-22U-S/M/L		

Compact Conical Abutment	
	STRAIGHT COMPACT CONICAL ABUTMENT DRIVERS
Handheld	I-AD
Handpiece	I-HAD
Wrench Insert	I-WI-A

Clinical benefits

The clinical benefits of the device are not dissimilar to that of dental implant therapy. These include improved chewing function, speech, aesthetics and patient psychological wellbeing. Through this procedure patients can expect to have their missing teeth replaced and/or crowns restored.

The clinical benefits specific to drills include:

- Effective cutting of bone within acceptable time limits.
- Precise preparation of the implant site to the dimensions of the intended implant.
- Can remain effective after multiple uses to the extent stipulated by the manufacturer.

The clinical benefits of handpiece devices include:

- Successful achievement of respective intended function (e.g., implant insertion, drill extension, screw removal)

Storage, cleaning and sterilisation

Limitations on reusable items

A direct value for reusable instruments cannot be given. Frequent processing may have minor effects on the instruments. The product life is normally determined by wear and damage during use, thus instruments if properly cared for and inspected after each use, can be reused many times. Maintain a checklist for these instruments recording the number of uses.

Prior to reprocessing the device, it should be thoroughly inspected and tested to determine its suitability for reuse.

NOTE: during use, handle drills and instruments with sterile tweezers to minimise contamination of the instrument tray and risk of damage to sterile surgical gloves.

Containment

As soon as practically possible, remove all visible residue after use (bone, blood or tissue), by immersing the instrument in cold water (dried soil may be difficult to remove).

Precleaning

Disassemble instruments from hand pieces and all connecting parts from instruments in order to clean soil from obstructed areas. Remove the PEEK bits from placement tools. Rinse with lukewarm water for 3 minutes and remove hardened debris with a soft nylon brush. Avoid mechanical damage during cleaning.

Manual cleaning or automated cleaning

Prepare an ultrasonic bath with suitable detergent (i.e. Steritech instrument cleaner - 5% dilution), sonicate for 20 minutes (alternative methods can be used if proven by the end user). Rinse with purified/sterile water.

NOTE: always follow the instructions for use of the manufacturers for cleaning agents and disinfectants.

Load devices into a thermodisinfectant. Run the cleaning and disinfection cycle, followed by the drying cycle.

Drying

Dry the instruments both inside and outside with filtered compressed air or single use, lint free wipes. Pack the instruments as quickly as possible after removal into the storage container. If additional drying is necessary, dry in a clean location.

Inspection

Do a visual inspection of the items to check for damages on instruments.

Packaging

Use the correct packaging material as indicated for steam sterilisation to ensure sterility is maintained. Double packaging is recommended. Where appropriate, the cleaned, disinfected and checked devices can be assembled and placed in instrument trays as applicable. Instrument trays can be double wrapped or placed in sterilisation bags.

Sterilisation

Southern Implants® recommends one of the following procedures to sterilise the restoration prior to use:

1. prevacuum sterilisation method: steam sterilise the abutments at 132°C (270°F) at 180 - 220 kPa for 4 minutes. Dry for at least 20 minutes in the chamber. Only an approved wrap or pouch for steam sterilisation must be used.
2. for users in the USA: prevacuum sterilisation method: wrapped, steam sterilise at 135°C (275°F) at 180 - 220 kPa for 3 minutes. Dry for 20 minutes in the chamber. Use a wrap or pouch that is cleared for the indicated steam sterilisation cycle.

NOTE: users in the USA must ensure that the steriliser, wrap or pouch, and all steriliser accessories are cleared by the FDA, for the intended sterilisation cycle.

Storage

Maintain packaging integrity to ensure sterility in storage. Packaging should be completely dry before storage to avoid corrosion and degradation of cutting edges.

Contraindications

Since the device is required or utilized as part of dental implant surgery, the contraindications for the device are not dissimilar to dental implant therapy as a whole.

The contraindications for dental implant therapy include:

- who are medically unfit for dental implant procedures
- where adequate numbers of implants could not be placed to achieve full functional support of the prosthesis,
- who are allergic or have hypersensitivity to pure titanium or titanium alloy (Ti-6Al-4V), gold, palladium, platinum, iridium or stainless steel.
- who are under the age of 18, have poor bone quality, blood disorders, infected implant site, vascular impairment, uncontrolled diabetes, drug or alcohol abuse, chronic high dose steroid therapy, anti-coagulant therapy, metabolic bone disease, radiotherapy treatment.

Warnings and precautions

THESE INSTRUCTIONS ARE NOT INTENDED AS A SUBSTITUTE FOR ADEQUATE TRAINING

- For the safe and effective use of dental implants, new technologies/systems and the device, it is strongly suggested that specialised training be undertaken. This training should include hands-on methods to gain competency on the proper technique, biomechanical requirements and radiographic evaluations required for the specific system.
- Improper technique can result in implant failure, damage to nerves/vessels and/or loss of supporting bone.
- Use of the device with incompatible or non-corresponding devices can result in device failure or poor performance.
- When handled intraorally, devices must be adequately secured against aspiration. Aspiration of products may lead to infection or physical injury.
- Electrosurgery should not be attempted around metal implants as they are conductive.
- The use of non-sterile items can lead to secondary infections of the tissue or transfer infectious diseases.
- Failure to perform appropriate cleaning, resterilisation and storage procedures as per the Instructions for Use document can result in device damage or secondary infection.
- Blunt drills may cause damage to the bone which could compromise osseointegration.
- For short implants, clinicians should closely monitor patients for any of the following conditions: peri implant bone loss, changes to implant's response to percussion, or radiographic changes in bone to implant contact along the implant's length. If the implant shows mobility or greater than 50% bone loss, the implant should be evaluated for possible removal. If the clinicians choose a short implant, then clinicians should consider a two-stage surgical approach, splinting a short implant to an additional implant, and placement of the widest possible fixture. Allow longer periods for osseointegration and avoid immediate loading.

NOTE: Training should be undertaken by both new and experienced implant users before using a new system or attempting to do a new treatment method.

It is important that thorough screening of prospective implant candidates is performed. This screening should include:

- a comprehensive medical and dental history.
- visual and radiological inspection to determine adequate bone dimensions, anatomical landmarks, occlusal conditions and periodontal health. Failure to accurately estimate lengths of drills relative to radiographic measurements can result in permanent injury to nerves or other vital structures.
- the evaluation and consideration of bruxism and unfavourable jaw relations.

Proper preoperative planning with a good team approach between well trained surgeons, restorative dentists and lab technicians is essential for successful implant treatment.

For successful implant treatment, it is important to:

- minimise the trauma to the host tissue as this increases the potential for successful osseointegration.
- be aware and avoid damage to vital anatomical structures such as nerves, veins and arteries. Injury to these structures may cause serious complications, including, injury to the eye, nerve damage and excessive bleeding. It is essential to protect the infraorbital nerve.
- identify actual measurements relative to the radiographic data as failing to do so could lead to complications.

Special care should be taken when treating patients with local or systemic risk factors that could affect the healing of the bone and soft tissue or otherwise increase the severity of the side effects, risk of complications and/or implant failure. These factors include:

- smoking.
- history of periodontal disease.
- poor oral hygiene.
- history of orofacial radiotherapy.**
- bruxism and unfavourable jaw relations.

** The potential for implant failure and other complications increases when implants are placed in irradiated bone as radiotherapy can result in progressive fibrosis of vessels and soft tissue, leading to diminished healing capacity.

The responsibility for proper patient selection, adequate training and experience in implant placement and the provision of the appropriate information required for informed consent rests with the practitioner.

Should the device not operate as intended, it must be reported to the manufacturer of the device. The contact information or the manufacturer of this device to report a change in performance is: sicomplaints@southernimplants.com

Side effects

Since the device is required or utilized as part of dental implant surgery, the side effects of the device are not dissimilar to those of dental implant therapy as a whole. Common side effects of dental implant therapy include pain, inflammation, phonetic difficulties and gingival inflammation. Other, less common side effects or complications to dental implant therapy include, but are not limited to: (1) wound dehiscence; (2) peri-implantitis; (3) transient weakness, numbness, and/or pain associated with mild nerve injury; (4) fat emboli formation; (5) marginal bone loss within acceptable limits; (6) allergic reaction(s) to implant and/or abutment material; (7) unspecified infection; (8) implant failure due to insufficient levels of osseointegration; (9) breakage of the implant and/or abutment; (10) loosening of the abutment screw and/or retaining screw; (11) complications requiring revision of the dental implant; (12) nerve damage resulting in permanent weakness, numbness, or pain; (13) loosening of the implant requiring revision surgery; (14) perforation of the maxillary sinus; (15) perforation of the labial and/or lingual plates and (16) bone loss possibly resulting in revision or removal of the implant.

Precaution: maintaining sterility protocol

Implants are packaged as follows:

1. An outer package consisting of a rigid, clear box which acts as protection for the inner package.
2. The inner package consisting of a blister pack (clear plastic-formed blister base with a TYVEK "peel-back" lid).
3. Within the inner package, there is a hollow tube which contains one implant suspended from a titanium ring, this ensures the implant never touches the inside of the plastic tube.
4. Labelling information is located on the surface of the peel-back lid and on the outside of the rigid box.

Care must be taken to maintain the sterility of the implant by proper opening of the packaging and handling of the implant.

1. Open the implant package in the non-sterile field, with non-sterile gloves, tear the address label to open the box.
2. With non-sterile gloves, remove the inner blister pack. Do not place the plastic box or blister pack-lid onto the sterile field. The contents of this inner package are sterile.
3. The sealed blister is to be opened by an assistant (with nonsterile gloves): remove the TYVEK lid and drop or place the sterile tube onto the sterile field, open the tube cap and attach the implant placement tool onto the implant and carefully remove from the sterile tube. Do not touch the sterile implant.

Other sterile components are packed in a peel pouch or blister base with a “peel-back” lid. Labelling information is located on the bottom half of the pouch, inside the packet or on the surface of the peel-back lid. Sterility is assured unless the pouch is damaged or opened. Non-sterile components are supplied clean but not sterile in a peel pouch or blister base with peelback lid. Labelling information is located on the bottom half of the pouch or on the surface of the peel-back lid.

Notice regarding serious incidents

Any serious incident that has occurred in relation with the device must be reported to the manufacturer of the device and the competent authority in the member state in which the user and/or patient is established.

The contact information for the manufacturer of this device to report a serious incident is as follows: sicomplaints@southernimplants.com.

Materials

Handpiece latch drivers:	Stainless Steel (DIN 1.4197)/SAE 316
Hand held drivers:	Stainless Steel (DIN 1.4197)/SAE 316
Idler:	Radel PPSU
Wrench drivers:	Stainless Steel (DIN 1.4197)/SAE 316

Disposal

Disposal of the device and its packaging: follow local regulations and environmental requirements, taking different contamination levels into account. When disposing of spent items, take care of sharp drills and instruments. Sufficient PPE must be used at all times.

Summary of Safety and Clinical Performance (SSCP)

As required by the European Medical Device Regulation (MDR; EU2017/745), a Summary of Safety and Clinical Performance (SSCP) is available for perusal with regard to Southern Implants® product ranges.

The relevant SSCP can be accessed at <https://ec.europa.eu/tools/eudamed>.

NOTE: the above website will be available upon the launch of the European Database on Medical Devices (EUDAMED).

Disclaimer of liability

This product is part of the Southern Implants® product range and should only be used with the associated original products and according to the recommendations as in the individual product catalogues. The user of this product has to study the development of the Southern Implants® product range and take full responsibility for the correct indications and use of this product. Southern Implants® does not assume liability for damage due to incorrect use. Please note that some Southern Implants® products may not be cleared or released for sale in all markets.

Basic UDI

Product	Basic-UDI Number
Basic-UDI for Reusable Instruments	6009544038769E
Basic-UDI for Drills and Handpiece Instruments	6009544038759C

Related literature and catalogues

CAT-2004 - Tri-Nex Implants Product Catalogue

CAT-2005 - IT Implants Product Catalogue

CAT-2010 - Osseointegrated Fixtures Catalogue

CAT-2020 - External Hex Implants Product Catalogue

CAT-2042 - Deep Conical Implants Product Catalogue

CAT-2043 - Internal Hex Implants Product Catalogue

CAT-2060 - PROVATA® Implants Product Catalogue

CAT-2069 - INVERTA® Implants Product Catalogue

CAT-2070 - Zygomatic Implants Product Catalogue

CAT-2093 - Single Platform Implants Product Catalogue

Symbols and warnings

											
Manufacturer: Southern Implants® 1 Albert Rd, P.O Box 605 IRENE, 0062, South Africa. Tel: +27 12 667 1046	CE mark 2797	Prescription device*	Sterilised using irradiation	Non-sterile	Use by date (mm-yy)	Do not reuse	Do not resterilise	Catalogue number	Batch code	Medical device	Authorised representative in the European Community
											
Authorised representative for Switzerland	Date of manufacture	Magnetic Resonance conditional	Magnetic Resonance safe	Single sterile barrier system with protective packaging inside	Single sterile barrier system	Consult instruction for use	Caution	Keep away from sunlight	Do not use if package is damaged		

* Prescription device: Rx only. Caution: Federal Law restricts this device to sale by or on the order of a licenced physician or dentist.

Canada licence exemption: Please note that not all products may have been licensed in accordance with Canadian law.

All rights reserved. Southern Implants®, the Southern Implants® logotype and all other trademarks used in this document are, if nothing else is stated or is evident from the context in a certain case, trademarks of Southern Implants®. Product images in this document are for illustration purposes only and do not necessarily represent the product accurately to scale. It is the responsibility of the clinician to inspect the symbols that appear on the packaging of the product in use.

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